


Connect with computer.

Wireless battery charge.


Jaundice Data Management Software by USB.

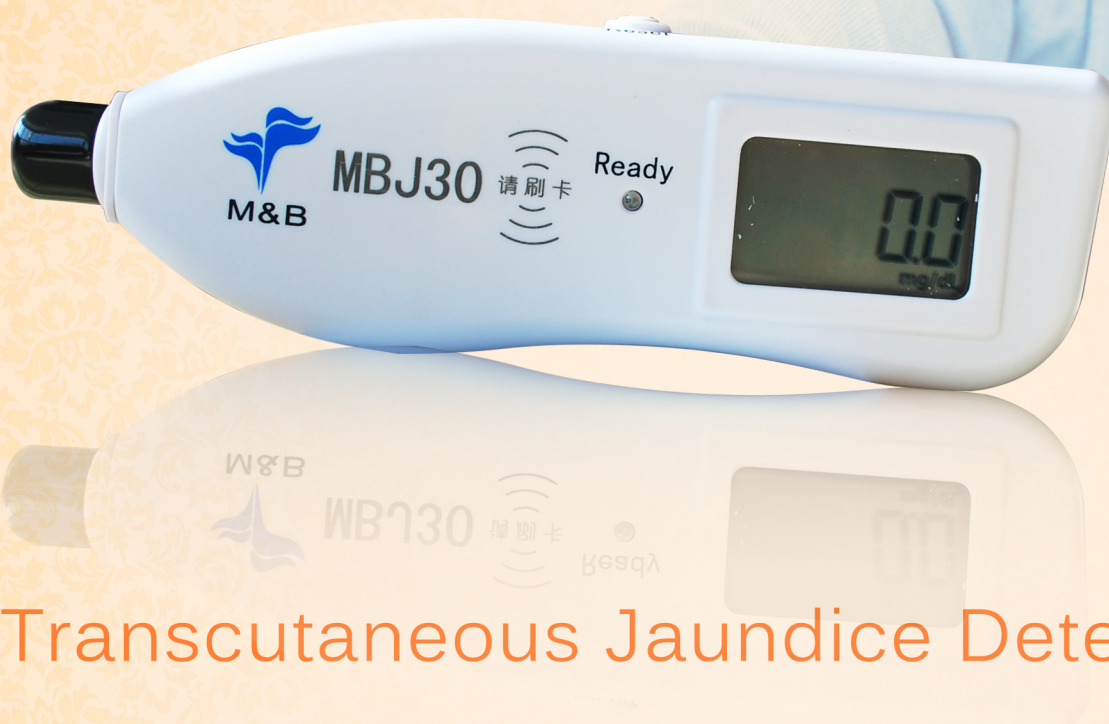

## Transcutaneous Jaundice Detector

- Wireless battery charge, Battery thermal protection
- LCD display for instant test reading
- Storage and memory function: 192 readings
- Patient ID card
- Average Times Setting from 1 to 9
- Upload the stored Jaundice Data to computer by USB



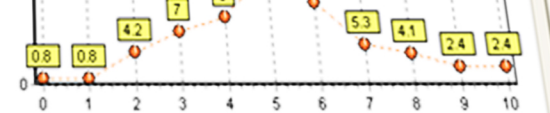

删除全部

显示全部


پرتوژن  
SHENASA PARTO GENE



MBJ30 is a handheld instrument used in the dynamic clinical examination of neonate jaundice. It tests transcutaneous bilirubin that is correlative with serum bilirubin concentration non-invasively and instantly through placing it on neonatal skin, which will prevent the neonates from any discomfort brought by blood drawing. Meanwhile it reduces test costs and increases the working efficiency of doctors and nurses.

## Specifications

- Display: LCD, 3 figures.
- Power: AA 1.2V×2 rechargeable battery.
- Indicator for ready: Green.
- Measurement range: 0.0mg/d L ~ 32.0 mg/d L.
- Measurement accuracy: Range is  $\pm 1.5$  mg/d L ( $\pm 25.5$   $\mu\text{mol/L}$ ).
- Charging preparation time: < 5 seconds.
- Record function: record 192 previous measurement values and review recorded data.
- Repeatability: < 10% .

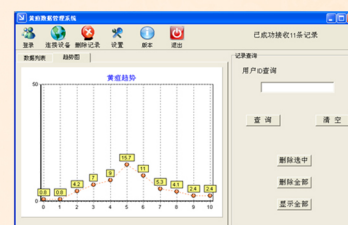

序号	用户ID	测量值	测量时间
176	499253132	0.8	2017-08-15 10:33:51
177	499253132	0.8	2017-08-15 10:34:00
178	499253132	4.2	2017-08-15 10:34:11
179	499253132	7	2017-08-15 10:34:21
180	499253132	9	2017-08-15 10:34:48
181	499253132	15.7	2017-08-15 10:52:22
182	499253132	11	2017-08-15 10:52:36
183	499253132	5.3	2017-08-15 10:52:49
184	499253132	4.1	2017-08-15 10:53:19
185	499253132	2.4	2017-08-15 10:53:38
186	499253132	2.4	2017-08-15 10:53:47

Beijing M&B Electronic Instruments Co., Ltd

Add:No.27,Yongwang Road, Beijing Bioengineering and Medicine Industry

Base, Huangcun Town, Daxing District, Beijing 102600 P.R, China

Connect with computer

Transcutaneous Jaundice Detector

